




## System Accessories

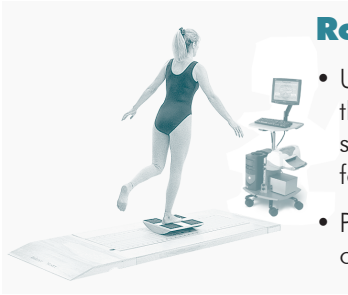


### Rocker board:

- Unstable surface to minimize the effectiveness of ankle strategy and right/left push-off forces for COG control.
- Provides unidirectional challenge of the feet/ankles bilaterally.

### Inversion/eversion wedges:

- 3° and 6° M/L
- Provides uneven surfaces for closed kinetic chain positioning.

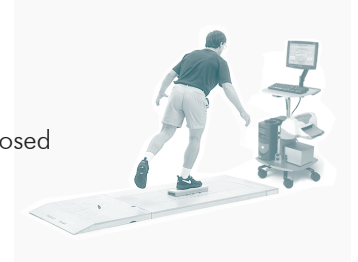

### Step-up blocks:

- 4 in. and 6 in. block
- Used in unilateral concentric/ eccentric strengthening exercises.

### Foam pad:

- 18 x 18 x 5 in (46 x 46 x 13 cm)
- 3.75 lbs. per cubic foot
- Unstable compliant surface to reduce somatosensory inputs for balance control during assessment and training.





### Leveling block:

- 2 in. block
- Used with wedges to level the pelvis and avoid leg length discrepancy effects.


### Heel/toe wedges:

- 6° and 12° A/P
- For stretching/strengthening of gastrocnemius and anterior tibialis muscles.
- Useful for patients with plantar flexion contractures.

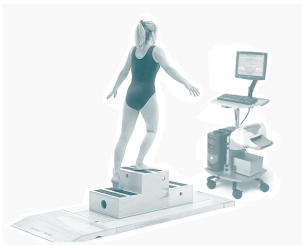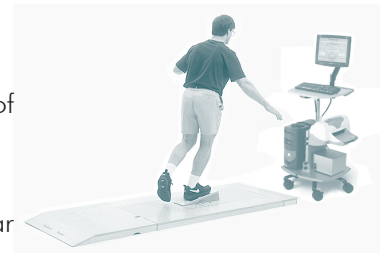

### Four piece step/stair set:

- 4 in., 8 in. and 12 in. step, and two-step chair
- Used in bilateral concentric/ eccentric strengthening exercises.

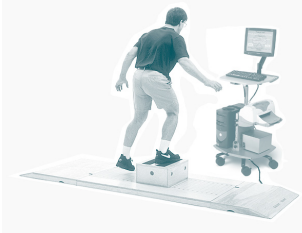

### Stability bar and harness:

- 46 in. x 46 in. x 86 in.
- Used with the static forceplate to provide a safer performance environment.





## System Accessories

